

Sterile disposable neoprene gloves, offering robust hand protection and fewer allergic reaction risks

- **Assured cleanliness:** DermaShield™ 73-711 sterile gloves are specially designed for compatibility with aseptic and Class 100 (ISO 5)/Grade A cleanroom environments
- **Reduced allergy risks:** These cleanroom gloves are also free from latex proteins and chemical accelerators, reducing wearer susceptibility to Type I or Type IV allergic reactions
- **Added durability:** Thanks to their unique neoprene formulation, these cleanroom safety gloves remain durable and resistant to punctures, while offering chemical splash protection
- **Enhanced features:** They also boast a beaded cuff with SUREFIT™ Technology, providing a more secure, stable fit and preventing cuff roll-down



KEY FEATURES & BENEFITS

- **Assured cleanliness:** Class 100 (ISO 5)/Grade A cleanroom compatibility
- **Latex-free, accelerator-free design:** Low allergy risks
- **Robust neoprene formulation:** Durability and chemical splash protection

Industries

- Biotechnology Manufacturing
- Medical Device Manufacturing
- Laboratory Research and Development
- Controlled and Critical Environments
- Production and Manufacturing

Recommended For

- Transferring liquids
- Clean room cleaning and preparing
- Sample taking and processing



TECHNICAL DATA SHEET


Product Information	
Material	Neoprene (Polychloroprene)
Color	Green
Shape	Anatomic with Curved Fingers
Cuff	Beaded with SUREFIT™ Technology
Manufacturing/QMS Audit Standards	ISO 9001:2015
Regulatory/Standard Compliance	AS/NZS 4179, ASTM D3577, ASTM D7160, CE 0493, EN 556-1:2001, EN ISO 21420:2020, EN ISO 374-1:2016, EN ISO 374-5:2016, EN 421:2010, ISO 10282, ISO 11193, ISO 13485, UKCA
Packaging	<p>1 pair per inner poly pack; 10 inner poly packs per inner polybag; 5 outer polybags per bag; 4 bags per master bag; 1 master bag of 200 pairs per carton/case</p> <p>More sustainable packaging: These sterile gloves are packed in recyclable plastic packaging* and delivered in recycled cardboard shipper cases.</p> <p>*Inner wrap, pouch, bag and liner are made from polyethylene (PE) based film. Always check your local recyclable status as these materials may not be considered suitable for recycling in your location.</p>
Storage	Keep away from direct sunlight; store in a dry place and keep in the original packaging. Keep away from ozone sources. If products are properly stored, as indicated, they won't lose their performances or change characteristics significantly. If products could be affected by ageing or storage, the expiry date is mentioned on the packaging materials.
Country of Origin	Sri Lanka
Available sizes	6, 6.5, 7, 7.5, 8, 8.5, 9
Powder Content	Powder-Free
External Glove Surface	Textured Fingertips
Internal Glove Surface	Polymer Coated with DERMASHIELD™ Technology
Sterilization Method	GAMMA irradiation (25 kGy)
Cleanroom Class	Class 100/ISO 5 & EU GMP Grade A
Shelf Life	5 Years
Tested for use with Chemotherapy Drugs	Yes, in accordance with ASTM D6978 (Not listed in the US FDA 510k)
Protein Level	N/A: contains no natural rubber latex
Anti-static	No

Physical Properties		Testing Method
Typical Length (mm/in)	300 / 12	EN 420/ASTM D3767
Freedom from Holes	0.65 AQL	EN 455-1/ASTM D3577
Typical Particle Count ≥0.5µm (counts / cm ²)	< 3500	IEST-RP-CC005.4
Target Single Wall Palm Thickness (mm/mil)	0.16 / 6.29	EN 420/ASTM D3767
Target Single Wall Finger Thickness (mm/mil)	0.20 / 7.87	EN 420/ASTM D3767
Target Single Wall Cuff Thickness (mm/mil)	0.15 / 5.90	EN 420/ASTM D3767
Ultimate tensile strength (MPa) Before Aging	17	ASTM D412-06a
Ultimate tensile strength (MPa) During Aging	12	ASTM D412-06a
Force at Break (N) Before Aging	9	EN 455-2
Force at Break (N) During Aging	9	EN 455-2

ORDERING INFORMATION

SIZE	6	6.5	7	7.5	8	8.5	9
REORDER NO.	73711060	73711065	73711070	73711075	73711080	73711085	73711090

PERFORMANCE STANDARDS AND REGULATORY COMPLIANCE

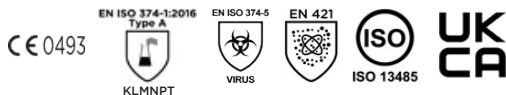

For additional information visit us at www.ansell.com, or call us at

Europe, Middle East & Africa Region

Ansell Healthcare Europe NV
T: +32 (0) 2 528 74 00
F: +32 (0) 2 528 74 01

Asia Pacific Region

Ansell Global Trading Center
T: +603 8310 6688
F: +603 8310 6699

North America Region

Ansell Healthcare Products LLC
US T: +1 800 800 0444
US F: +1 800 800 0445
CA T: +1-800-363-8340

Latin America & Caribbean Region

Ansell Commercial Mexico S.A. de C.V.
T: +52 442 248 1544 / 248 3133

Australia

Ansell Limited
T: +61 1800 337 041
F: +61 1800 803 578

UK

Ansell Nitritex
T: +44 1638 663338
F: +44 1638 668890

Technology



Cuff Roll-down Technology



Ansell, ® and ™ are trademarks owned by Ansell Limited or one of its affiliates. US Patented and US and non-US Patents Pending: www.ansell.com/patentmarking © 2025 Ansell Limited. All Rights Reserved.

Neither this document nor any other statement made herein by or on behalf of Ansell should be construed as a warranty of merchantability or that any Ansell product is fit for a particular purpose. Ansell assumes no responsibility for the suitability or adequacy of an end user's selection of gloves for a specific application.